
**JOHNSON®**

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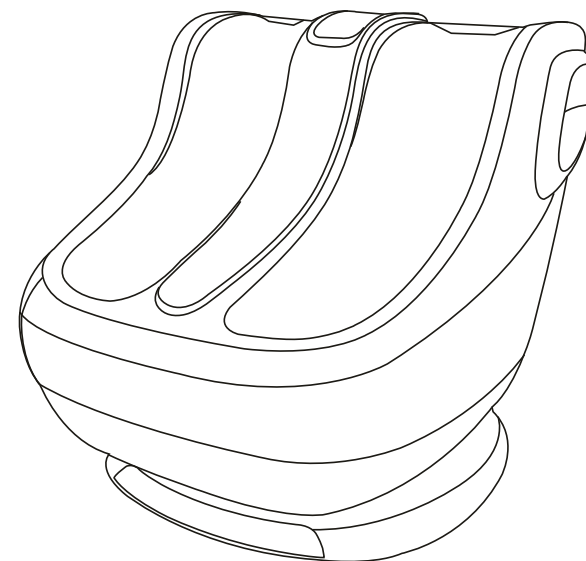Johnson Health Tech., Co., Ltd.  
[www.johnsonfitness.com](http://www.johnsonfitness.com)

**足部按摩器**  
**Leg & Foot massager**

**RT1889**  
**使用手冊**  
USER MANUAL

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The name and content of hazardous substances in this product are as follows

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Statement of environmental protection use period

Use this product under the described condition of user manual, the period that will not leak harmful substances or elements, not to cause environment pollution or personal and property damage is 10years.

## ► TROUBLE SHOOTING

NO.	Symptom	Possible Cause	Solution
1.	There is a clattering noise in the movement.	This is perfectly normal caused by the working of product.	It is normal and no need for treatment.
2.	The Foot Massager can not be started up.	The power plug is loosened.	Insert the plug into socket firmly.
		The power switch is not turned on.	Turn on the power switch.
3.	Stop during use	Default working time is up.	Reuse in 15 minutes later

If there are some unusual symptoms after having checked the above items, please contact your local dealer to ask for an authorized service engineer.

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## ► CLEANING AND MAINTENANCE

### Maintenance

1. Cleaning
  - Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
  - Use a soft cloth to wipe the surface of the product and controller.
  - Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of massage chair, and then use a piece of dry and soft cloth to wipe it.

### ⚠ Caution

- Do not pull the power plug by damp hands.
- Do not use water washing product, do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the massage chair, otherwise, the color of upholstery will fade or the massage chair will crack or be eroded.
- Do not iron the upholstery.
- When product is wet, do not detach the cloth cover, let the product dry naturally.

2. Maintenance
  - Make sure the power supply is turned off and the plug has been pulled out from the socket.
  - Sweep off the dust. Keep the massage chair clean and away from damp places.
  - If the massage chair is not to be used for a long period of time, please cover the massage chair with a dustproof cover(selfprovided).

### ⚠ Caution

- Do not store the product under direct sunlight or in places with high temperature. This may cause the color of the upholstery to fade.

### Kindly reminding

Please check if there is any below symptoms when the product is in daily use:

- The burned and scorched flavor.
- Touch the wire, sometimes power is on, sometimes power is off.
- The power cord heats.
- Other abnormalities.


### ⚠ Caution

- If any of the above mentioned abnormalities exist, please stop use to avoid accident happens.
- Turn off the power supply and pull the plug out, contact the Caution local dealer or contact us.
- Do not disassemble the product and repair.

## ► USE METHOD

### Preparation before use

- Please put the product on the flat floor.
- Please keep this product one meter away from audio-visual equipment such as TV set, radio, etc to prevent any signal interference.
- Please plug to rated power supply(100-240~).
- Turn on the power switch.
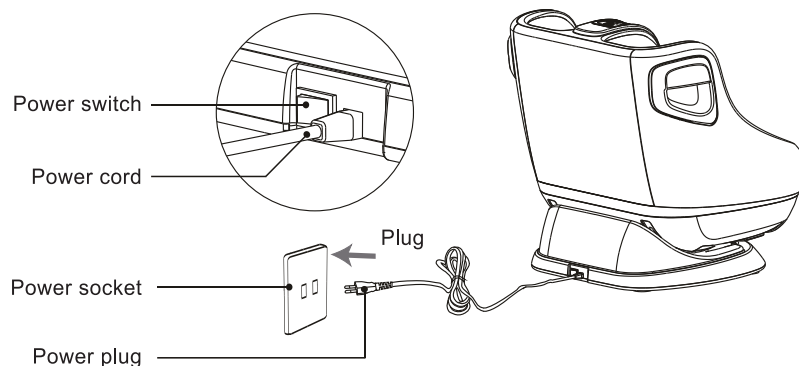
### Caution

- Before plug, please check and make sure the power cord is undamaged.
- Do not use the wet hand to pull out the power cord in order to avoid any electric shock.

- Please choose a chair with desired height and be seated, then put your feet into the product (see right picture).
- Please turn on the power button “”, start the auto massage program.
- During using, if want to adjust the angle of the product, please press the "angle adjustment button" (inside the right handle). Press this button up can adjust the angle the products. Three grades, max angle 22 degrees. (please see the picture on page 3)



### After use

- Turn off the power button on the control panel.
- Remove the plug from the power supply socket.

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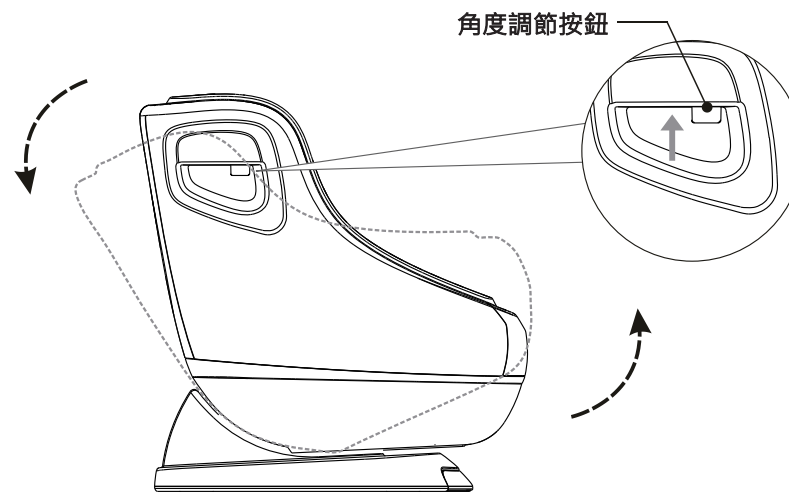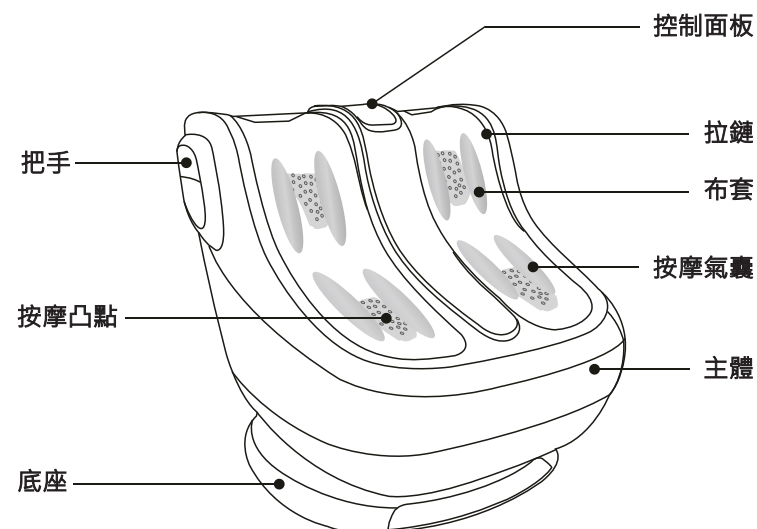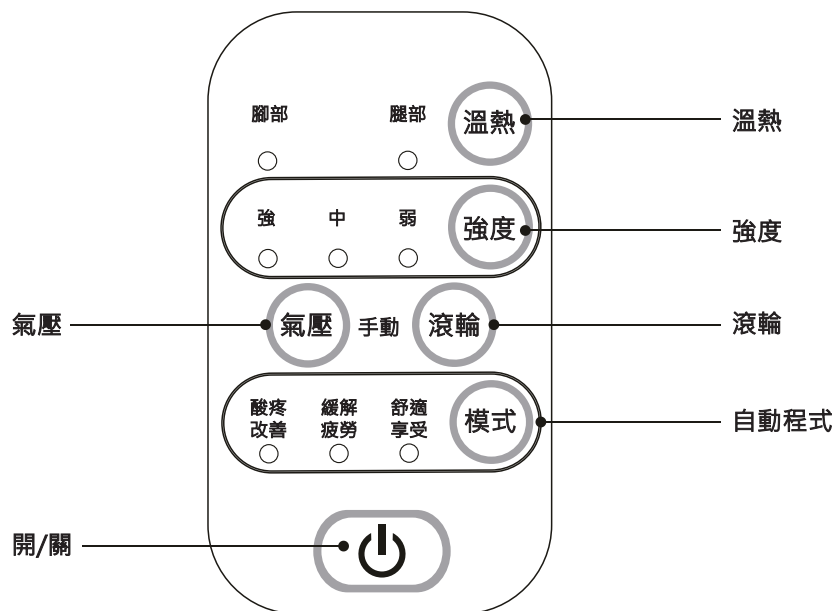
### Functions of control panel

#### Heat button

Press this button, cycle through to turn on the heat function for the calf, sole and turn off the heat function for the calf and sole. (Default setting: turn on the heat function for the calf and sole)

#### Intensity button

Press this button, adjust the intensity of scrape and air pressure. Press this button, cycle through to choose the strong, middle, weak grade. (The corresponding indicator will be lightened or extinguished).

#### Air pressure button

The air pressure is a manual mode, can turn on or turn off the air pressure massage. Press this button, start the air pressure massage function, no scrape. Indicator light flashes slow, the intensity of air pressure massage become slow. Indicator flashes fast, the intensity of air pressure massage become fast.

#### Scrape button

Scrape button is a manual mode. Can turn on or turn off this function. Press this button is scrape function, no air pressure massage function. Indicator light flashes slowly, the intensity of scrape become slowly. Indicator light flashes fast, the intensity of scrape become fast. (Default setting: turn on this function)

#### Auto mode button

Press this button, cycle through to select three modes, reduce pain relieve fatigue, comfortable enjoyment. (Default setting: turn on the comfortable enjoyment)

#### ON/OFF


After connecting the power, press this button to start the massage function of product. The corresponding indicator will be lightened. Press this button again, turn off the massage function of the product, the corresponding indicator will be extinguished. (Default setting time: 15 minutes)

### ⚠ Caution

- ▶ When power failure happens during massage, please turn off the power and pull out of the plug from output socket.
- ▶ Please do not let your hand or fingers close to and inserted into the gap of products for preventing any injuries caused.
- ▶ Please do not let children and pets play the products for preventing any accident happened.

## ▶ NAME AND FUNCTION OF EACH PART

### Control Panel
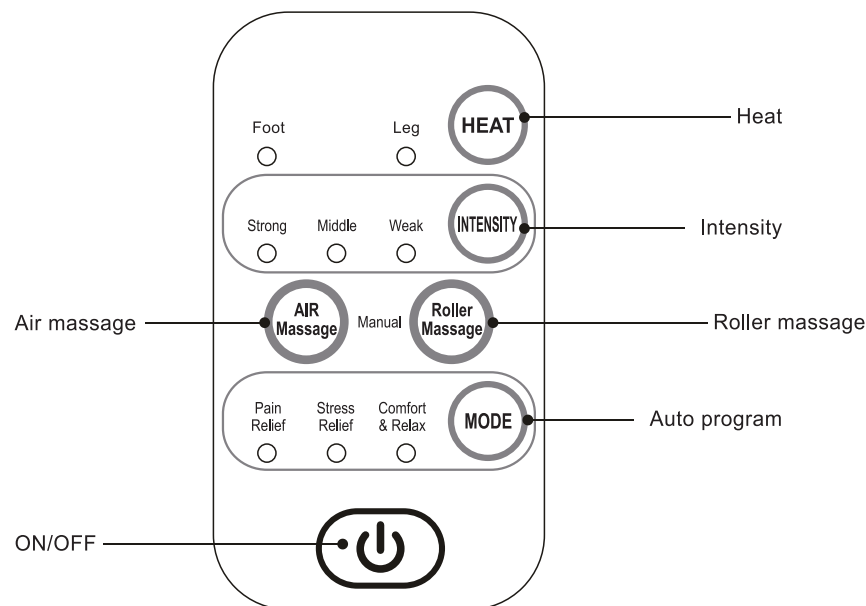
### ⚠ Caution

- ▶ Please don't put water or other liquid on the control panel.
- ▶ Please don't press the button too long time at a time to avoid the program confused.
- ▶ Please don't use the product by standing posture or put any heavy object on the product to avoid any broken caused.
- ▶ Please don't stand or seat on the product.
- ▶ When moving the product, please don't pull the power cord or drag the product directly to avoid any damaged caused.

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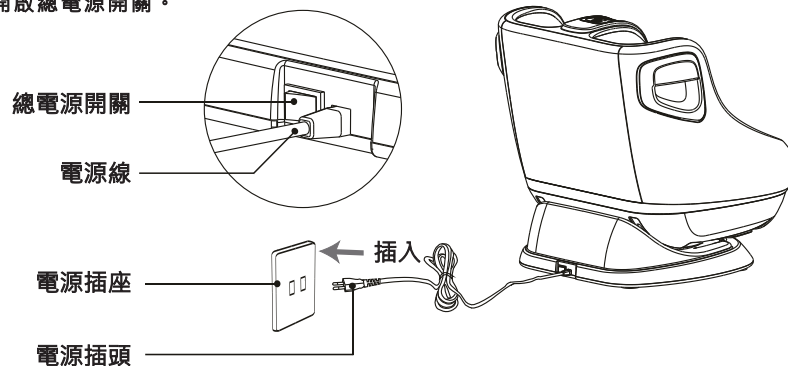
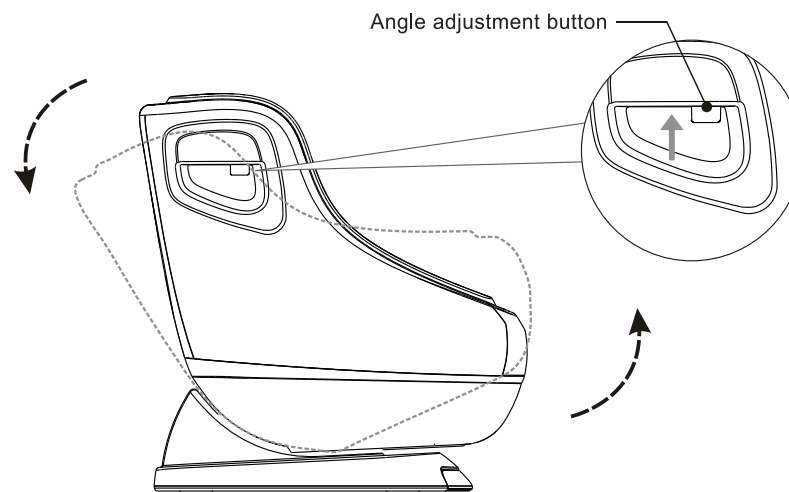
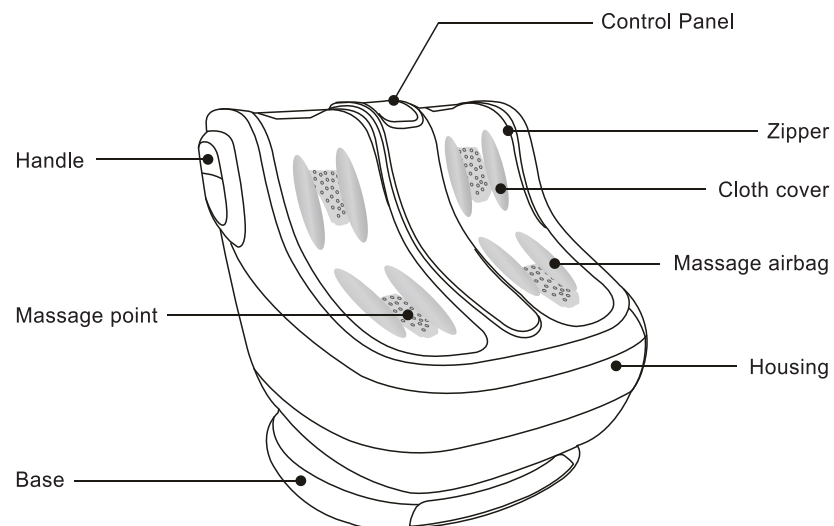
## ► SAFETY PRECAUTIONS

### ⚠ CAUTION

COMPULSARY	Do not continuously use this product over 15 minutes.
	Do not put heavy objects on the product.
	The default working time of the product is 15 minutes.
	Please do not cover the product with blanket or pillow when it is used, thus to avoid user from being burn or causing fire disaster.
	Please do not stand on the product to avoid any failure.
	Do not operate the controller with wet hands.
	Do not put the product outdoors.
	It is forbidden to dip the whole product in water or other liquid when cleaning.
	If the product is moved from a low temperature place to a warmer place, it is advised to use it one hour later to avoid function failure caused from vapor condensing.
	Do not use this product together with other therapeutic equipments.
	People who are insensitive to heat must be accompanied by normal person when using this product.
	Do not use this product for medical treatment.
	Please stop using this product immediately when you find any abnormal situations (such as smoke or peculiar smell) during use.
	Please use the power supply with rated voltage only.
	Please check the product before using, if damage found, please contact the dealer for repairing immediately.
PROHIBITED	Do not use this product in damp environment to avoid electric shock.
	During use, if you feel abnormal, please stop using this product immediately, and then consult your doctor.
	After use, please cut off power supply.
	Do not use the product within half an hour after meal.
	Do not sleep during the product working.
	Do not use the product after drunk or any abnormality of your body.
	Do not wear watch, bangle and not put cell phone, keys and other valuables in your jacket when use the product.
	Do not use petrol, caustic agents and spay pesticide for cleaning.
	Do not put the product near combustible.
	Do not smoke when using the product to avoid fire occurrence.

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## ► SAFETY PRECAUTIONS

### Safety Precautions

- Before using, please read this "SAFETY PRECAUTIONS" carefully to ensure proper operation.
- The following "SAFETY PRECAUTIONS" is of great help to your correctly use of this product.
- To clearly indicate the degrees of danger and damage, the classes of accidents caused by incorrect use are classified into "WARNING" and "CAUTION", instructions marked as followings should be strictly followed.

	<b>WARNING</b>	Improper operation may cause injury to users.
	<b>CAUTION</b>	Improper operation may cause users injury or product damage.
		Forbidden. Specific instruction will be indicated in the diagram drawings. (e.g, the left drawings means disassembly is not allowed.)
		Compulsory Specific instruction will be indicated in the diagram drawings. (e.g, the left drawings means to disconnect the power plug.)

### WARNING

 <b>COMPULSARY</b>	We recommend that those with any of the following conditions consult a doctor before using.
	<div> 1. Those who are using electronic medical instrument imbedded in the body, such as cardiac pacemaker etc.  2. Those who are being treated by doctors, especially those who feel uncomfortable.  3. Those who are malignant tumor patients, heart disease patients, acute patients. </div> <div> 4. Women who are pregnant or menstruating.  5. Those with osteoporosis or fracture spine.  6. Those with skin disease or those whose skins are injured.  7. Those whose body temperature is over 38 degrees(Febrile Phase) </div>
 <b>PROHIBITED</b>	Supervision is compulsory when the product is used by intellectual disabilities children.
	Do not allow children play with product.
	Do not more than one person to use this product at a time.
	Do not damage, process, excessively bend, pull hard, twist or knot power cord.
	Do not use the product when the power cord is damaged.
 <b>NO DISASSEMBLY</b>	Do not allow the minors to use the product.
	Anyone else, except maintenance personnel specified by our company, should not disassemble, maintain, or remodel this product. If any problem, please contact local dealer immediately.

Thank you for purchasing Johnson RT 1889 Leg & Foot massager. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start. This manual should be kept on hand for easy reference on safe and efficient operation of the Massager.

## ▶ CONTENTS

SAFETY PRECAUTIONS .....	12
NAME AND FUNCTION OF EACH PART .....	14
External Structure .....	14
Control Panel .....	15
Functions of Control Panel .....	16
USE METHOD .....	17
Preparation before Use .....	17
Operation Method .....	17
After Use .....	17
CLEANING AND MAINTENANCE .....	18
Maintenance .....	18
TROUBLE SHOOTING .....	19
TECHNICAL SPECIFICATIONS .....	20

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