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J Pro X Muscle Massage Gun Heat + Cool Therapy

JOHNSON®

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Customer Care Line at 0800-017360
www.johnsonfitness.com.tw

Customer Care Line at 1300 88 0558
www.johnsonfitness.com.my

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V. Recycling

Once the service life of the product ends, please do not discard it with the general household wastes. You should place the product to the recycle bin appointed by the government. It is also acceptable to send the product back to the after-sales service center of the original factory, and the personnel in the service center will dispose it in an environmentally sound manner.

1. Maintenance, cleaning, storage and transportation

01. Clean the device with a slightly wet towel and wipe dry with soft cloth.
02. Please store the device in the packing box, if it will not be use for a long time.
03. Please charge the device once a month, to prevent it from being affected from moist or humid condition.
04. The massage device should be stored in a place that have good ventilation and free from corrosive gas, with a distance of 250mm from the ground.
05. Do not exceed 5 layers when stacking the massage devices.
06. Please clean the massage device before packing & storage if it is not use for a long time.
07. The product is applicable to the transport mode including vehicle, train, ship and air transportation.
08. To prevent violent vibrations, collisions, rolls and falls, the product package should be placed upwards, handled with care and prevented contact from water.

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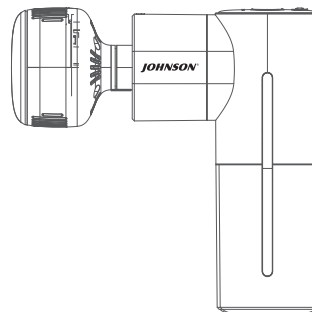
**注意：本手冊之英文說明內容並不適用於在台灣所購買之型號。
The explanations in the English version of this manual are not applicable to the models purchased in Taiwan.**

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TROUBLESHOOT	
PROBLEM	TROUBLESHOOT
1.The massage head slips off and drops.	Check whether the massage head is installed in place (see the manual for specific operations)
2. Device could not be turned on	Check whether the battery is full (please charge for more than 5 hours if it is not use for a long time)
3. Device could not be charged	Check whether the charging port is plugged in place correctly
4. Device suddenly stops after continuous running for 10 minutes	Restart the device by pressing the ON - OFF switch.
5.The knocking function stops after the massage head is installed	Restart the device by pressing the ON - OFF switch.

IV. FAQs

If the device fail to on:

1. Environment Temperature > 40°C

If the device operates in an environment of >40°C, the device may trigger overheat protection and will not start.

Solution:

Keep the device in use environment <40°C for 2 hours and then re-start it.

2. Battery power is at 0:

The device battery will produce a over-discharge protection if the device is not charged for a long time.

Solution:

Fully charge the device using the cable provided and start as normal.

If the devise fail to charge:

1. Device Temperature too High

In case of high temperature detected in the main body of the device, it will not charge even if a charger is accessed.

Solution:

Keep the device in an environment under 40°C to cool down for 2 hours and then charge as normal.

2. Excessive temperature of the charger.

If the use environment is above 40°C, the device may trigger the temperature overheating protection. It will not charge even if a charger is accessed.

Solution:

Keep the device in an environment under 40°C to cool down for 2 hours and then charge as normal.

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
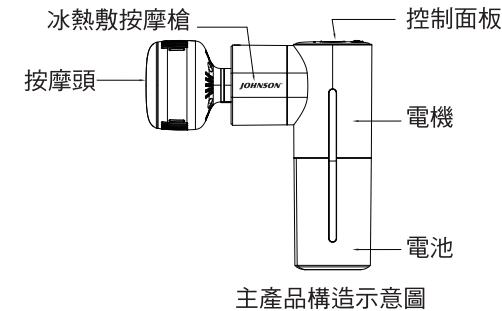
- 03. Active charging is indicated when battery power LED lights flash regularly.
- 04. Full charge is indicated when battery power LED lights stop flashing and remain fully illuminated.
- 05. Battery can be recharged at anytime and at any battery level.
- 06. Recommend time to recharge when battery level power is lower than 25%.
- 07. It is not recommended to completely run down the battery level to 25% before charge.
- 08. The average run-time is 3 hours depending on speed level and applied pressure during use.

III. Basic parameters

Product name: J PRO X Muscle Massage Gun	Product model: Johnson C66-M
Product net weight: 620g	Input voltage: 5V DC, 2A
Power: 25W (maximum 35W)	Battery specification: 11.1V 2400mAh
Strength level: 3 levels available for adjustment	Charging time: About 5 hours
Service environment: 0-40°C	Host size: 140×60×150mm
Battery life: 2-3 hours for knocking 2-3 hours for cooling 3-4 hours for heating	Self-timing: Knocking for 15 minutes Cooling for 15 minutes Heating for 15 minutes
Temperature value: Heating:35°C for Level-1, 40°C for Level-2, 45°C for Level-3 Cooling:20°C for Level-1, 15°C for Level-2, 10°C for Level-3 There will be a plus or minus deviation of 5°C in the above values, depending on the environment and service conditions.	


Massage head user manual:

1. Plug and unplug the massage head to assemble or change.
2. Fully insert the massage head in the piston notch.

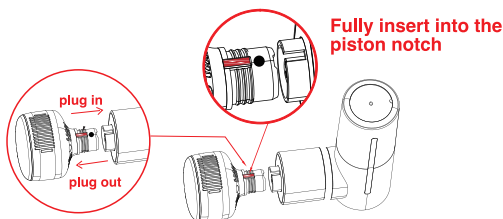
Charge

The charging port is located in the top end of the product handle lever.

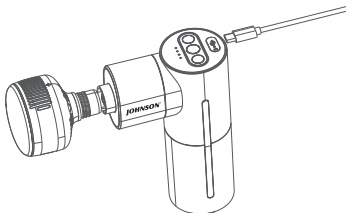01. Fully charge the battery up to 5 hours before first use.
02. This device utilize a Type C cable for charging
(provided together in the package)

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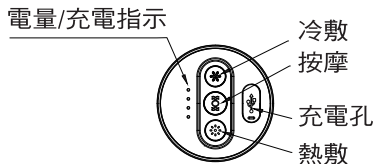4. In power OFF mode: long press the COOL Therapy switch for 2 seconds to turn on the device to Level 1 temperature of 20°C (blue light), press twice for Level 2 temperature of 15°C (purple light), press three times for Level 3 temperature of 10 °C (red light). Long press the COOL Therapy switch again for 2 seconds to turn the device off.
5. In the power ON mode: press any function switch to activate the desired function mode.
6. The device will enter standby mode after the OFF switch is press. The battery indicator will continue display for 5 seconds. Press any function switch during this period to restart the device. After 5 seconds of inactivity, the product will automatically shut down.
7. For the HEAT+COOL therapy, the HEAT+COOL massage head must be installed to activate the function. Other massage heads are not applicable.
8. Vibration, HEAT and COOL therapy mode to be used separately.
9. Vibration, HEAT and COOL therapy mode automatically timed for 15 minutes per session.

Massage head assembly

HEAT+COOL head

Round head

Bullet head


Flat head

U Shaped head

Note: The yellow solid line in the figure is the product alignment line, and it is black on device.

II. Operating manual

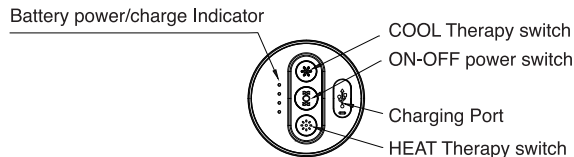Control Panel Diagram

1. In power OFF mode: press any of the function key to show current battery power status.
2. In power OFF mode: long press the ON-OFF power switch for 2 seconds to turn the device on to Level 1 intensity (blue light), press twice for Level 2 intensity (purple light), press three times for Level 3 intensity (red light). Long press the ON-OFF power switch again for 2 seconds to turn the device off.
3. In power OFF mode: long press the HEAT Therapy switch for 2 seconds to turn the device on to Level 1 temperature of 35°C (blue light), press twice for Level 2 temperature of 40°C (purple light), press three times for Level 3 temperature of 45°C (red light). Long press the HEAT Therapy switch again for 2 seconds to turn the device off.

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
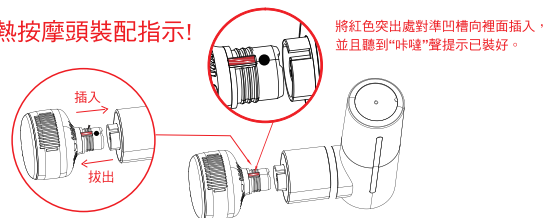
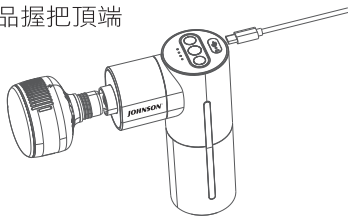
I. Specification

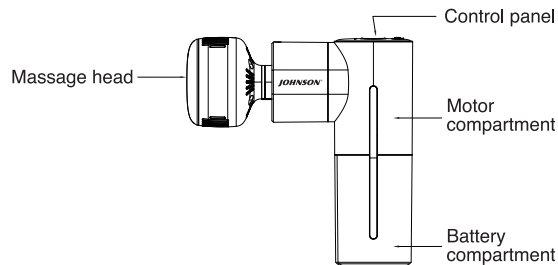Massage head

HEAT+COOL head
Mainly used for relaxing your whole body, shaping and skincare

Bullet head
Mainly used for deep tissue relaxation & precise massage

Round head
Suitable for spine, large muscle groups & thigh muscle massage

Flat head
Mainly used for large contact area, bulk muscle groups and shaping

U Shaped head
Mainly used for the neck, cervical spines and heel tendon

Diagram: Do not touch this part to avoid pinching when the device is started without a massage head.

Diagram: Do not touch this part to avoid pinching when the device is installed with a massage head.

2. Disclaimer

Please do not use the massage gun or any vibration device without first obtaining approval from your doctor if any of the following apply.

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown the integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

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Important Information

To reduce dangers including fault, electric shock, injury, fire, death, and device damage, please read all instructions, warnings and reminders of the product in this operation manual carefully before using the device.

1. Safety Instruction

01. For adult use only.
02. Do not use the device on the wet or dirty skin surface. In this case, you can use it wearing your clothes by gently pressing and moving on the skin, each position for about 60 seconds.
03. Do not immerse the product in the water or allow water enter into the device ventilation hole.
04. Do not place any objects into the ventilation hole of the device.
05. Do not put the device into fire or close to any fire source.
06. Do not intentionally drop the device.
07. Do not tamper or modify the device in anyway.
08. Never leave the device operating or charging unattended.
09. Keep fingers, hair or any body parts away from the shaft and back side of the device as pinching might occur.
10. Due to the heat-producing surface of the device, please use the device carefully.
11. Do not use the device during charging mode.


Thanks for purchasing our products. Please carefully read this operation manual before use to ensure correct operation. Please keep this operation manual well for future reference.

Brief Introduction and Purpose of Massage Gun

J Pro X Muscle Massage Gun is a device designed for deeply relaxing muscular fascia and suitable for a post-workout massage. It can alleviate stiffness and ache of muscles, relieve muscle fatigue, and accelerate blood circulation to improve the soft tissue health of the body.

Purpose

Suitable for physical massage and skincare before/after exercise and fitness.

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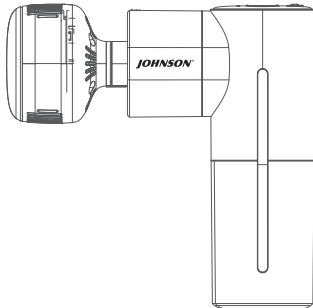Declaration of the Presence Conditions of the Restricted Substances Marking

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